

Software Market Release Notes

eviFluor Duo Fluorometer Instrument Firmware and API Software



Instrument manufacturer:

Hombrechtikon Systems Engineering AG (HSE)

Document number:

DMR-DOU-21200611-04962-007_A2

Release date:

26.01.2026

Table of Contents

1	Document change history	3
2	About This document	3
3	Instrument Firmware	4
3.1	Version History	4
3.2	What's New in Version 0.8.0	4
3.3	Update Recommendation & Impact	4
4	API Software Interface C / CLI.....	5
4.1	Version History	5
4.2	What's New in Version 0.6.0	5
4.3	Update Recommendation & Impact	5
5	API Software Interface C#	5
5.1	Version History	5
5.2	What's New in Version 0.12.0.....	5
5.3	Update Recommendation & Impact	6
6	Compatibility	6
6.1	Instrument Firmware / Instrument Hardware	6
6.2	Instrument Firmware / API Software Interface.....	6
7	Instructions for Instrument Firmware update.....	7
7.1	Required tools.....	7
7.2	Instrument Firmware update procedure	7

1 Document change history

Document version	Release date	Change
007	26.01.2026	• API Software Interface C / CLI Version 0.8.0 • API Software Interface C# 0.10.0 • Addition of Chapter "Compatibility" • Addition of Chapter "Instructions for instrument firmware update"
006	22.08.2025	• First customer facing version

2 About This document

The eviFluor Duo Fluorometer is a compact instrument designed for integration into an automatic liquid handling system.

This document contains the revision history and changes for the eviFluor Duo Fluorometer software packages and is intended to give a detailed overview for liquid handler software developers as well as for liquid handler sales and service representatives.

- Instrument firmware: The instrument firmware runs on eviFluor Duo Fluorometer instruments and is responsible for the actual optical measurement.
- API Software interface: The API Software interface can be integrated into the liquid handler software and is responsible for obtaining the fluorescence values and the calculation of the nucleic acid concentration. The API is available in the following programming languages:
 - C / Command line interface (CLI)
 - C#

The API software can be accessed via github: <https://github.com/hseag/evifluor>

- Upon request, the API can be made available in Python as well

Incidents and support requests shall be reported to the HSE eviFamily Service Desk:
<https://hseag.atlassian.net/servicedesk/customer/portals>

3 Instrument Firmware

3.1 Version History

Release version	Release date	Changes / Updates
0.8.0	11.06. 25	Initial release for end customer usage: - Internal logging functionality added
0.6.0	28.03.25	Autogain function for LED power setting improved for performance optimization
0.4.1	03.03.25	Support of updated version of hardware PCBA
0.4.0	25.02.25	Initial release (internal)

3.2 What's New in Version 0.8.0

- Release date: 11.06.2025
- Resolved issues: Not applicable (initial release)
- Improvements / New features: Internal logging functionality only
- Known issues: None

3.3 Update Recommendation & Impact

- Update recommendation: Not applicable (initial release)
- Impact: Not applicable

4 API Software Interface C / CLI

4.1 Version History

Release version	Release date	Changes / Updates
0.6.0	26.01.26	Addition of Run class, in order to facilitate implementation/ integration on liquid handlers
0.4.0	25.02.25	Initial release

4.2 What's New in Version 0.6.0

- Release Date: 26.01.26
- Resolved Issues: None
- Improvements / New features: Addition of Run class, to facilitate implementation/ integration on liquid handlers
- Known issues: None

4.3 Update Recommendation & Impact

- Update recommendation: Update not required by reseller (HSE internal functionality only)
- Impact: Impact: No impact on reseller integration nor instrument performance

5 API Software Interface C#

5.1 Version History

Release version	Release date	Changes / Updates
0.12.0	26.01.26	- Functionality added to report progress of firmware update to eviManager Software - Addition of Run class, in order to facilitate implementation/ integration on liquid handlers - C# implementation made runnable on Linux for HSE unit tests or other future needs
0.10.0	21.08.25	Additional check method added to verification class to check if concentration result is negative.
0.8.0	11.06.25	Kit-based concentration fit added -> Parameter added to differentiate between High-Sensitivity (HS) and Broad Range (BR) Assay kits. Currently only HS supported.
0.6.0	14.03.25	Added "Get Available Devices" method
0.4.0	25.02.25	Initial release

5.2 What's New in Version 0.12.0

- Release Date: 26.01.2026
- Resolved Issues: None
- Improvements / New features:
 - The eviFluor Duo C# API is used in the new eviManager Software. In order to be able to display a progress bar on the eviManager Software during a firmware update, additional functionality was added to eviFluor Duo API C#
 - Addition of Run class, to facilitate implementation/ integration on liquid handlers

- C# implementation made runnable on Linux for HSE unit tests or other future needs
- Known issues: None

5.3 Update Recommendation & Impact

- Update recommendation: **An update to version 0.12.0 by resellers is not required**, as there is no new functionality which impacts field usage. The added functionality is for HSE internal purposes only.
- Impact: Impact: No impact on reseller integration nor instrument performance

6 Compatibility

The following subchapters describe the compatibility between the firmware versions, API versions and the instrument hardware.

6.1 Instrument Firmware / Instrument Hardware

The current instrument firmware is compatible with all released instrument hardware versions.

6.2 Instrument Firmware / API Software Interface

All API versions for C and C# are compatible with all instrument firmware versions.

7 Instructions for Instrument Firmware update


7.1 Required tools

In order to update the instrument firmware, the following is needed:


- eviManager Software
- PC / Laptop
- USB connection cable: USB Type-A to USB Type-C

7.2 Instrument Firmware update procedure


1. Download the latest version of the **eviManager Software** from GitHub:
<https://github.com/hseag/evimanager>
2. Download the latest firmware version from GitHub and save it to your PC.
<https://github.com/hseag/evifluor/releases>
- see respective release, file with extension .srec
3. Connect the eviFluor Duo instrument to your PC using a USB Type-A to USB Type-C connection cable. (USB Type-C socket on bottom of eviFluor Duo instrument)

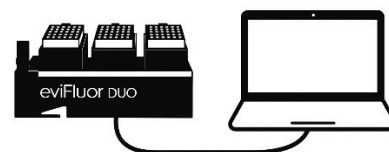


4. Open the **eviManager Software** on the PC: The software will automatically detect the connected instrument and display the serial number and current firmware version in the "Information" tab.

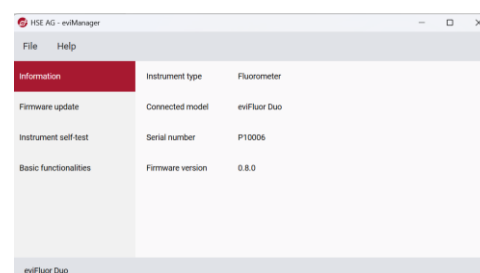


5. Click on the tab "Firmware update"

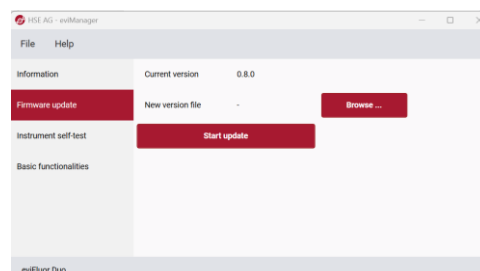


- Click "Browse" and select the downloaded firmware file

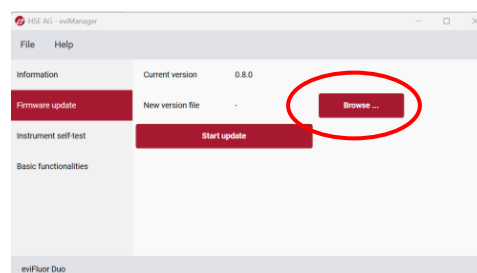

- Click on "Start update"

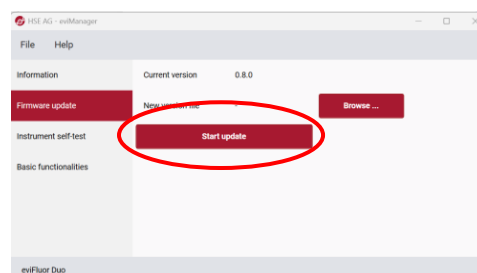

- The update will automatically begin and indicate the progress on the button.

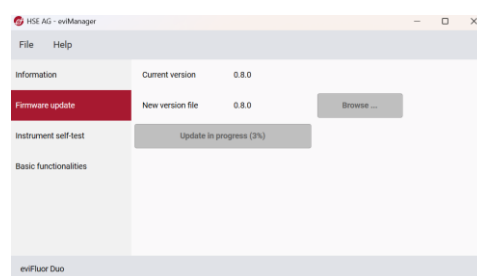

- Once the update is complete, an on-screen message will confirm successful completion of the update.

